**bébé·confort®**

**Thermomètre auriculaire**  
**Ear thermometer**  
**Termometro auricolare**  
**Oorthermometer**  
**Termómetro de oído**  
**Termómetro auricular**  
**Лобный термометр**

**3106203400**

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
## EN - Ear Thermometer

Please be sure to read all warnings, instructions and notations before using this product. Please keep these instructions as they contain important information.

### WARNING:

- **NEVER** allow a child to play with this thermometer. It is not a toy. It is a sensitive electronic device and may be damaged if not used properly under parental supervision.
- **ALWAYS** supervise your child when using thermometer.
- Avoid sudden impact and dust.
- **NEVER** allow your child to use it alone or without close parental supervision.
- **DO NOT** allow child to walk or run during temperature taking. Probe could present a hazard.
- **NEVER** immerse this device in water or other liquids(not waterproof).

### 1 - Know Your Unit
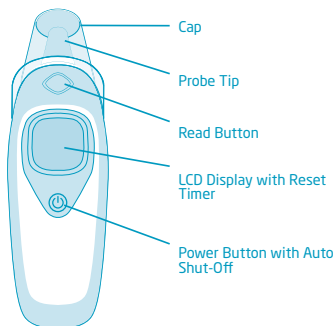
### How Does it Work?

The Easy Read Ear Thermometer obtains a reading in 1 second by measuring the infrared heat given off by the eardrum and surrounding tissues.

### Temperature From The Ear

Clinical research has shown that the ear is an ideal site for taking body temperature. The eardrum shares the same blood supply with the hypothalamus, the part of the brain that controls body temperature. Therefore, the ear is an accurate indicator of core body temperature.

### 2 - Important Information

Read these instructions carefully to ensure accurate temperature readings. Temperature readings are affected by many factors, including exercise, drinking hot or cold beverages prior to taking a reading, and technique.

- High or prolonged fever requires medical attention, especially for young children. Contact your physician. Fever chart included on package.
- Clean probe tip before and after each use with a damp cloth or alcohol swab.
- Wait 60 seconds for thermometer to reset between readings.
- Keep in mind that the thermometer will automatically shutoff after one minute in ready mode.
- Do not expose the thermometer to extreme temperatures, humidity, direct sunlight or shock. Store the thermometer in a dry place at room temperature.
- Damaged device is not reliable. Do not use the device if you think it is damaged or if anything appears unusual.
- Do not use device in the electromagnetic environment. Using the unit in the immediate vicinity of mobile phones, microwave appliances or other devices with strong electromagnetic fields may result in impaired functioning.
- Batteries and electronic instruments must be disposed of in accordance with the locally applicable regulations, and not as domestic waste.
- Do not use this device close to strong electromagnetic fields such as mobile telephones or radio installations. Keep a distance from such devices when using this unit.

### 3 - Positioning Thermometer

For an infant, it is best to have the child lying flat with his head sideways so the ear is facing upwards. For an older child or adult, it is best to stand behind and slightly to the side of the child. It is recommended that three temperatures in the same ear be taken 60 seconds apart and the highest one used as the reading.

- Proper placement within the ear canal is essential for an accurate reading. If the thermometer is not inserted properly and firmly into the ear canal, a lower reading might be obtained. Remember to always use the same ear. The temperature between the right and the left ear may differ.

### 4 - Operating Thermometer

1- **Turn On:** Press and release the Power Button. All segments are displayed for 1 second (Figure 1). The last temperature measured appears with backlight for 2 seconds (Figure 2).

2- Thermometer is ready to take temperature when it displays 3 lines and flashing "F/C" and diamond on screen (Figure 3). Read button will also flash green.

3- Position ear for thermometer:

- **Under 1 Year:** Gently pull ear straight back

- **1 Year to Adult:** Gently pull ear up and back (See "Positioning Thermometer")

4- Insert thermometer until the probe tip seals the ear canal.

5- Press the Read Button. After 1 second thermometer will beep once and display temperature.

**FEVER ALERT NOTE:** If patient has a fever, thermometer will sound multiple beeps.

6- Reset timer will count up to 60 seconds to indicate next time thermometer will be ready to take a reading. Once reset timer is complete, thermometer will return to ready mode and will auto shut off after another 60 seconds (Figure 4).

NOTE: You can halt this process at any time by pressing the power button to turn the thermometer off, and then turning it back on again.

7- Auto-shut off will occur after 1 minute, or press Power Button to turn off.

All segments

Fig.1

Previous Temp.

Fig.2

Ready to Read

Fig.3

Reset Timer

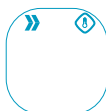
Fig.4

### 5 - Memory Function

The thermometer's memory stores the last reading.

It will appear after thermometer runs its built-in test at startup. This feature allows you to compare your most recent measurement to the next one.

### 6 - Care and Cleaning

- Very important! To ensure accurate readings, clean the probe and tip with an alcohol swab or cotton swab before and after each measurement.
- Use a soft, dry cloth to clean thermometer body. Never use abrasive cleaning agents, thinners or benzene for cleaning.
- Do not scratch the surface of the probe lens or the display.
- Do not submerge any part of the product in water.
- Make sure the thermometer has been at room temperature for at least 30 minutes before use.

## 7 - To Replace Battery (CR2032 or equiv.)

When the "low battery" icon appears (Figure 1), the battery is nearly exhausted, but a few measurements may still be taken. When "low battery" icon begins flashing, the battery is exhausted and needs replacing.

1 - Use any small, flat object (like a coin) to unlock the battery door (Figure 2).

2 - Press release tab with finger to remove battery (Figure 3). If battery does not easily remove, use wooden toothpick to help lift it from battery compartment.

3 - Insert a new 3-volt battery (CR2032 or equivalent), correctly positioned with the positive "+" side facing up (Figure 4) and gently press into place.

4 - Place the battery cover on the thermometer and secure it using a coin.

NOTE: Remove battery if stored for a long period of time.




CAUTION: Replace only with 3V "CR2032" or direct equivalent lithium battery. Use of other types of batteries may present a risk of fire or rupture.

### WARNING:

- Keep battery out of child's reach. Swallowing battery could be fatal. Battery should not be charged or placed into extreme heat as it may explode.
- Used batteries should be disposed of properly. Do not dispose of batteries in fire, recharge, put in backwards, disassemble, mix with used or other battery types - they may explode or leak and cause injury.

**END OF LIFE:** At the end of its service life, please dispose of the battery and thermometer at appropriate collection sites in accordance with currently applicable regulations.
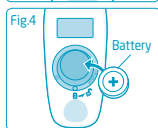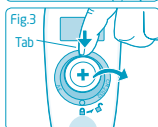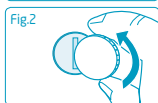
## 8 - Changing from Fahrenheit to Celsius

To switch the display between °F and °C: While the thermometer is off, press and hold the Power Button for about 3 seconds until you see the °F or °C. Release button.

## 9 - Troubleshooting

Problem	What to Do
Display reads <b>HI</b>	<ul style="list-style-type: none"><li>• Measured temperature is above 42.2°C. Turn off thermometer. Wait one minute and then turn on. Listen for two beeps before placing probe in ear. Re-take temperature</li></ul>
Display reads <b>Err</b>	<ul style="list-style-type: none"><li>• Thermometer has been in a hot or cold environment and must be allowed to reach room temperature. Wait at least an hour before attempting again to use this thermometer</li></ul>

Temp reading doesn't seem accurate

- Thermometer's probe tip may be dirty. Clean with alcohol swab and retake temperature.
- Ear canal may have excessive wax build-up. Clean ears per your pediatrician's advice.
- Install new battery and retry function before calling Consumer Relations on electronics.

### Technical Specifications

Measurement Range: 32°C to 42.2°C.

Measurement Accuracy: Technical: +/- 0.2°C.

Operation Environment: 10°C to 40°C, 15-95 % relative maximum humidit.

Storage Environment: -25°C to 55°C .15-95 % relative maximum humidity.

a.Device standard: Device corresponds to the requirements of the standard for infrared thermometers:

ASTM E 1965 EN12470-5 IEC 60601-1 IEC 60601-1-2 IEC60601-1-11

b.Electromagnetic Compatibility: Device fulfills the stipulations of the standard IEC 60601-1-2

c.The stipulations of EU-Directive 93/42/EEC for Medical Devices Class IIa have been fulfilled.

- The battery life. Battery lifetime:approx.500 measurements
- Expected service life:5 years
- IEC 60601-1-11

**IP22** IP22: Protected against solid foreign objects of 12.5 mm diameter and greater. Protected against vertically falling water drops when the device is tilted up to 15°.

### Warranty

This product is warranted for a two-year period from the purchase date. However, this warranty does not apply to damages, in particular, resulting from non-compliance with this guide, incorrect manipulations, shocks or inappropriate installation of the battery. It does not apply either to the battery or the package.

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Microlife Corporation, 9F, 431, RuiGuang Road, NeiHu, Taipei 11492, Taiwan, R.O.C.
MANDATAIRE / MANDATORY: MICROLIFE AG
**DOREL FRANCE S.A.S.** - Z.I. - 9 bd du Poitou, BP 905, 49309 Cholet Cedex, FRANCE
**DOREL SWITZERLAND S.A.** - Chemin de la Colice 4, 1023 Crissier, SWITZERLAND / SUISSE